

■ LDP1000 SERIES

The smaller, the better

■ PLASMADETEK-E

New era of detection

■ PLASMADETEK 3

The new "HOT" products

LD MAG

2016 EDITION



MULTIDETEK 2

ENHANCED DESIGN AND NEW CAPABILITIES

Where **innovation** leads to **success**

LDetek

CHROMALYTIC +61(0)3 9762 2034
ECHnology Pty Ltd

Website NEW : www.chromalytic.net.au E-mail : info@chromtech.net.au Tel: 03 9762 2034 . . . in AUSTRALIA

Australian Distributors
Importers & Manufacturers
www.chromtech.net.au

1 (of 16) 2016Mar

LDETEK COMPANY PROFILE

LDetek is a Canadian based company specialized in manufacturing, developing and integrating online gas analyzer, gas chromatograph system and many other accessories. Our mission is to provide reliable and performing systems with the best technology on the market.

With its commitment to continuous improvement, LDetek thinking and members always push the technology limit further. It brings innovative products and solutions to a market in a severe need of new and better quality materials. The wide range of products, applications and patents developed, make LDetek product portfolio attractive to many different spheres of gas analysis activities. With a constant objective of offering quality products, LDetek always makes all efforts to get rigorous quality control and all necessary levels of approbation and/or certification.

Our products and services are provided by a team of specialists with a strong experience involved in gas analysis business.

With a well established network of offices and partners in over 35 countries, a complete and proactive worldwide support is provided.

CONTACT INFORMATION

Head office (Thetford Mines, Canada)

Sales and general request: info@ldetek.com

Technical support: support@ldetek.com

Phone : + 1 418 755 1319

Fax : + 1 418 755 1329

Chinese office

rli@ldetek.com

Phone: + 86-10-87269199

Fax: + 86-10-87249394

Mobile: + 86 13911632686



Partners Listing

Argentina, Buenos Aires	Nigeria, Lagos
Australia, Melbourne	Norway, Oslo
Australia, Peakhurst	Pakistan, Karachi
Brazil, Rio De Janeiro	Philippines, Laguna, Binan
Brazil, Sao Paulo	Poland, Wroclaw
Chile, Santiago	Romania, Bucharest
Colombia, Bogota	Russia, St-Petersburg
Egypt, Cairo	Singapore
Finland, Helsinki	South Africa, Edenvale
France, Paris	Spain, Madrid
Germany, Frankfurt	Sweden, Stockholm
Greece, Athens	Taiwan, Taipei
Hungary, Budapest	Thailand, Bangkok
India, Mumbai	Turkey, Istanbul
Indonesia, Jakarta	UAE, Abu Dhabi
Iran, Tehran	Ukraine, Kiev
Israel, Acre	United Kingdom, Oxon, Wantage
Italy, Cesano Maderno	USA, Colorado, Denver
Japan, Tokyo	USA, Oregon, Dayton
Korea, Seoul	USA, Pennsylvania, King of Prussia USA,
Malaysia, Selangor	Texas, Houston
Mexico, Monterrey	Việt Nam, Ho Chi Minh City
Netherlands, Lelystad	



Where innovation leads to success

www.ldetek.com

THE NEW HOT PRODUCT!

PLASMADETEK 3

The PlasmaDetek 3 (patent pending) is a new detector with unique capabilities not available on the market yet. It is the only commercialized heated plasma emission detector that exists for gas chromatograph. By keeping the same principle and features of the PlasmaDetek 2 (spectroscopic emission with a micro-plasma), this new design offers the following new capabilities:

- **Heating with up to 200 Celsius**
- **Removable/changeable optical filters for a more flexible detector**
- **Changeable cell**
- **Possibility to connect a spectrometer fiber optic direct to the cell for specific lab or research project**



Patent pending

The PlasmaDetek 3 uses the plasma controller of the PlasmaDetek 2 which makes it compatible with any previous installation.

Clarity from DataApex will become compatible with the PlasmaDetek3, avoiding the use of the analog signal and simplify his use and installation.



NEW ERA



OF LIGHT HYDROCARBON MEASUREMENT

Looking for a **SAFE, SENSITIVE, LOW OPERATION COST** and **MAINTENANCE FREE** system?

The patent pending **PlasmaDetek-E** is the solution



SAFE

No fuel (H_2) and no related safety accessories



SENSITIVE

< 1 ppb lowest detection possible



LOW OPERATING COST


Only N_2 carrier gas supply



MAINTENANCE FREE

No periodic detector cleaning

TYPICAL INSTALLATION FOR LIGHT HYDROCARBON MEASUREMENT

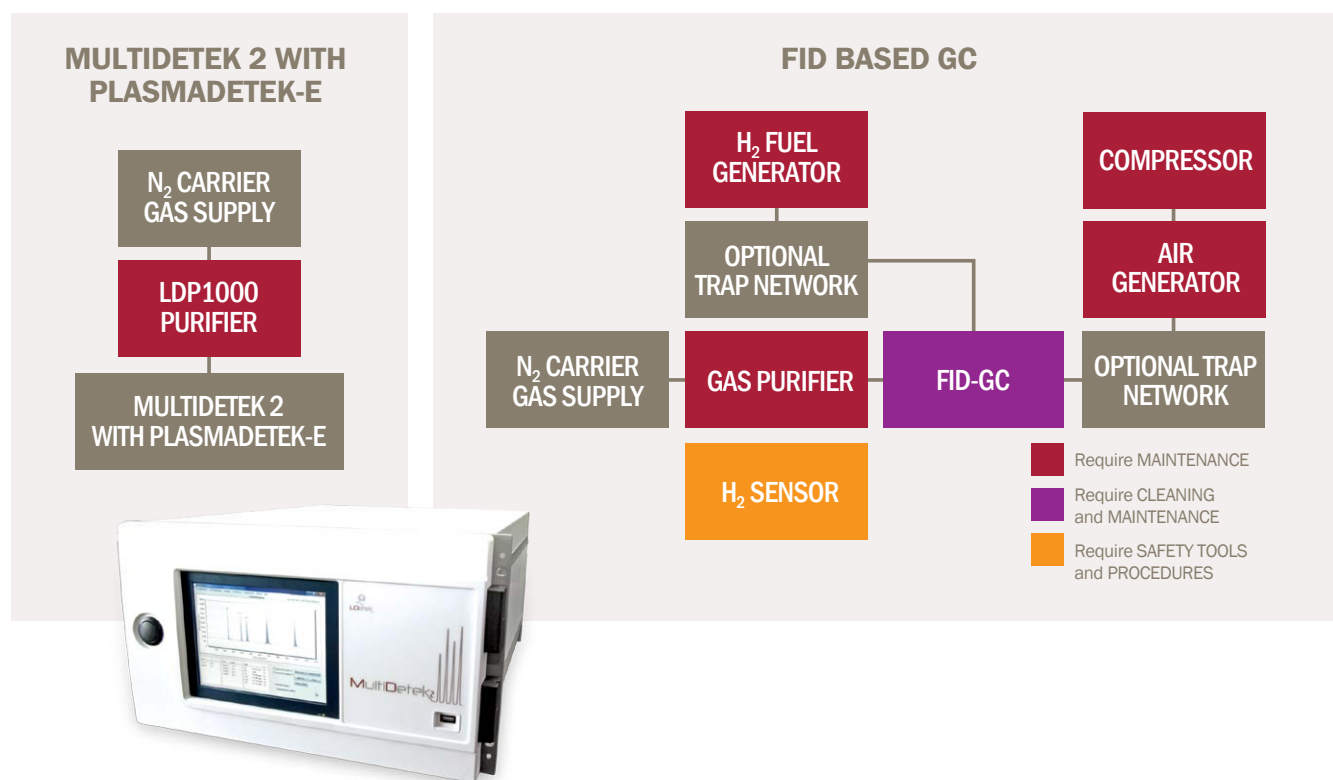

INSTALLATION COST* COMPARISON

Parts	MultiDetek 2 with PlasmaDetek-E	FID-GC
H ₂ generator	N/A	\$7200
Zero air generator	N/A	\$2125
Air compressor	N/A	\$1200
H ₂ safety accessories	N/A	\$1500
2 year maintenance cost	\$2000	\$5000
Total cost*	\$2000	\$17 025

* costs are approximate and may vary for each system


N/A: not applicable

Please consult Application Note LD14-01 on LDetek web site for more technical details.




PARTNER

Antelia integrated PlasmaDetek series solutions in many GCs over the last years and is proud of its strong expertise for developing high performance applications. Antelia works with Perkin Elmer lab GC platform to develop their solutions. On top of that, Antelia also integrates PlasmaDetek on the compact GAS GC. Depending on the applications, Antelia uses the GAS compact series or the Perkin Elmer desktop GC.

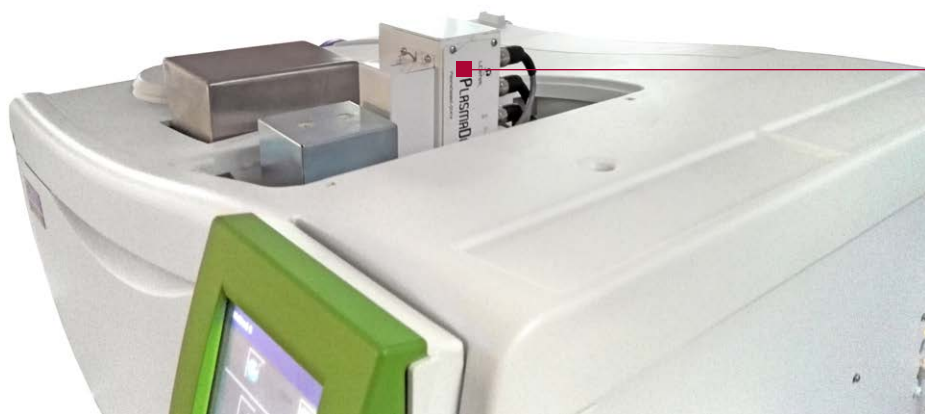

PlasmaDetek2(PED) mounted on the Clarus Perkin Elmer GC

PlasmaDetek1(PED) mounted on the Clarus Perkin Elmer GC

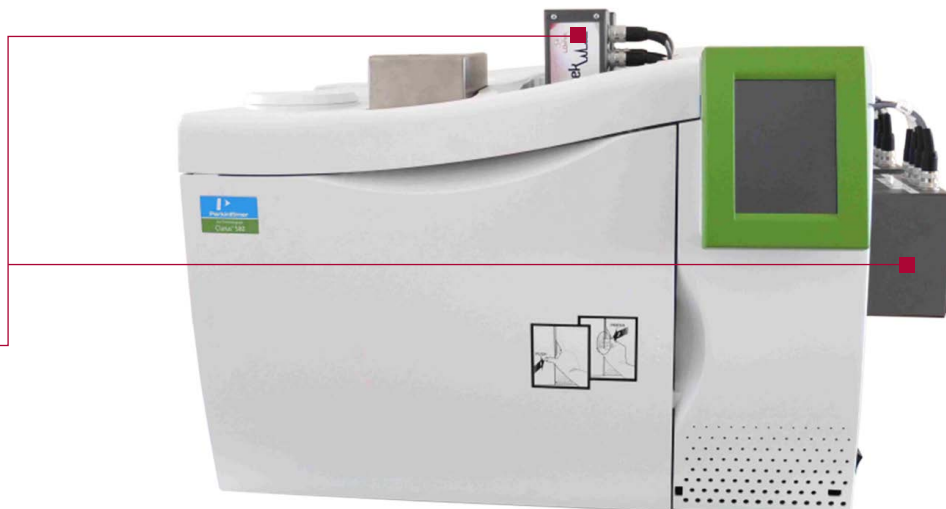

ANTELIA: 9 Chemin de la Brocardière, 69570 DARDILLY • contact@antelia.fr • Port +33 6 77 01 51 25 • Tel +33 4 72 17 82 10 • Fax +33 4 78 66 07 03

Where innovation leads to success

www.idetek.com

APPLICATION SUMMARY

As part of a research project for measuring greenhouse gases in France, more precisely for the N₂O and CO₂ for multiple soil sample analysis having a wide range of concentrations, the use of the PlasmaDetek2 (PED) detector and the Agilent 7890 GC combined with the Headspace 111 auto sampler system have been used in this proposed solution.

CHALLENGE OF THE APPLICATION

The samples from different locations come in 40ml vials. More than 4000 samples have to be analyzed per campaign. The samples contain CO₂ concentrations ranging from 500ppm up to 80% and N₂O concentration from 300ppb up to 5000ppm. The system of detection must be able to cover each sample within 15 minute cycle time. Low and high concentrations are analyzed with the same system.


SOLUTION

System implementation comprises a Headspace setter 111 positions coupled to an Agilent 7890 gas chromatograph (GC) equipped with purged valves and connected to two types of detectors, one TCD for high levels of CO₂ in particular, and a PED (PlasmaDetek2 from LDetek) for traces of N₂O. The advantage of this notorious PED detector with respect to ECD detector, frequently encountered in this application, is the absence of radioactive source subject to a heavy administrative procedure and staff empowerment. The presence of these two sensors in the same instrument, allows analysis and quantification of high levels of CO₂ and very low levels of N₂O.


SYSTEM PICTURE



RESULTS

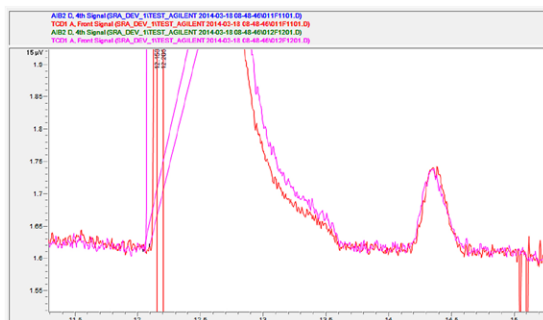

Chromatogram showing 80% CO₂ and 40ppm N₂O using TCD channel. This channel gives the ability to measure high concentration CO₂.

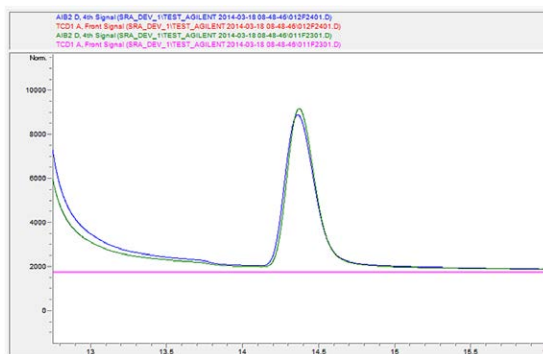

Chromatogram showing 40ppm N₂O using PED (PlasmaDetek2) channel. This channel is used for measuring low concentration N₂O.

SRAINSTRUMENTS
150 rue des Sources
69280 Marcy l'Etoile
France

T : 04.78.44.29.47
F : 04.78.44.29.62
info@sra-instruments.com
www.sra-instruments.com

SA à Directoire et Conseil de surveillance au capital de 150.000 €
RCS Lyon B 342 068 731
APE 4669B
SIRET: 342 068 731 00039
Code TVA FR 40342068731



Agilent Technologies

Premier Solution Partner

KOVALENT

DANI

A DANI PARTNER

PLASMADETEK2 (PED) MOUNTED ON DANI MASTER GC

APPLICATION SUMMARY

At Microbial Analytics Sweden AB we have investigated gases dissolved in water for over two decades. The water is collected from deep boreholes in solid rock and the gases are dissolved in the water under very high pressure. The background is that some of these gases might act as fuel for microorganisms; other gases indicate ongoing microbiological processes. It is important to investigate these mechanisms because they might endanger long-term storage of for instance spent nuclear fuel. Also, in a shorter time span other structures underground might be heavily corroded by such mechanisms.

SOLUTION

We have developed our own sampler design to collect and transfer the water pressurized to the laboratory. After extraction, we analyze the dissolved gases using several different gas chromatographs. The ranges of gases are the usual permanent gases, though present at unusual ratios and concentrations, as well as lighter hydrocarbons.

As a research institute it is important for us to keep the quality at a very high level. Because of this we have decided that all analytical results should be verified by multiple analyses on different gas chromatographs using alternative columns as well as alternative principles of detection.

Last year we decided to replace an older gas chromatograph that was set up with TCD detector as well as an FID detector and a methanizer. It was brought to our attention by the Swedish agent for PlasmaDetek2, Kovalent AB, that the plasma detector could be a useful alternative for the replacement instrument. Because of this we settled for the PlasmaDetek2 instead of the FID/methanizer setup on one channel in the new chromatograph, keeping TCD detection on the other channel.

For us this was a very fortunate decision since the PlasmaDetek2 has qualities difficult to obtain in other ways. Not only is it very easy to quantify carbon monoxide and

dioxide in the same run as the hydrocarbons without the need for a methanizer that has to be protected by switching valves. It is also possible to detect inert permanent gases in the same run. This is very much in line with our quality goals, to detect the analytes using different principles. Other instruments in our laboratory are still equipped with for instance FID detection and the results are compared, showing very good correlation.

Another positive effect was the ability to tune the sensitivity for different gases. In our samples we sometimes have to measure hydrogen at low concentration in the presence of comparatively high levels of neon. Chromatographic separation is possible, but not always optimal. Using the PlasmaDetek2, the sensitivity for hydrogen is much higher than for neon and the quantification gets a lot easier.

The PlasmaDetek2 (PED) detector has been a very positive experience for us!

Where innovation leads to success

www.ldetek.com

CHROMALYTIC +61(0)3 9762 2034
ECHnology Pty Ltd

Website NEW : www.chromalytic.net.au E-mail : info@chromtech.net.au Tel: 03 9762 2034 . . . in AUSTRALIA

Australian Distributors
Importers & Manufacturers
www.chromtech.net.au


RESULTS

The chromatogram below shows an example of analysis of a sample gas containing 0,10% of CO-CO₂-C₂H₂-C₂H₄-C₂H₆ and C₃ in balance nitrogen. The CH₄ is also measured in this application and the integration window appears at minute 5, but it is not present in the sample gas used for this analysis.


On the chromatogram, we can see that the CO peak is integrated in a slight drift that comes from the balance Nitrogen peak that elutes just before. The selectivity of the Plasmadetek2 (PED) for the CO measurement makes it suitable to measure it, even if the column used cannot separate N₂ from CO perfectly.

The C₂'s are also shown in a slight drift that comes from the temperature programming of the column. The peaks can be separated and integrated with success at any concentration.

The configuration of the system uses a single Carboxen 1010 column for the separation. The column temperature was programmed to allow the peaks having a late elution time to come earlier. The injection volume used was 100 micro liters.

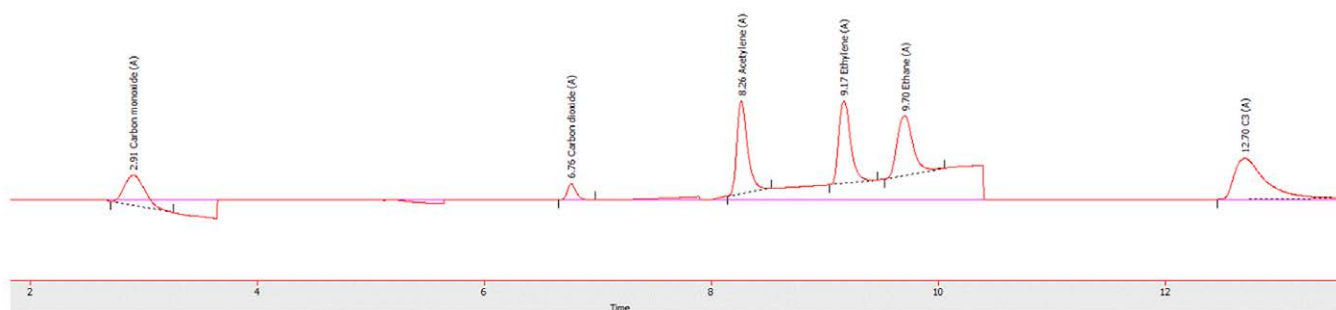

SYSTEM PICTURE

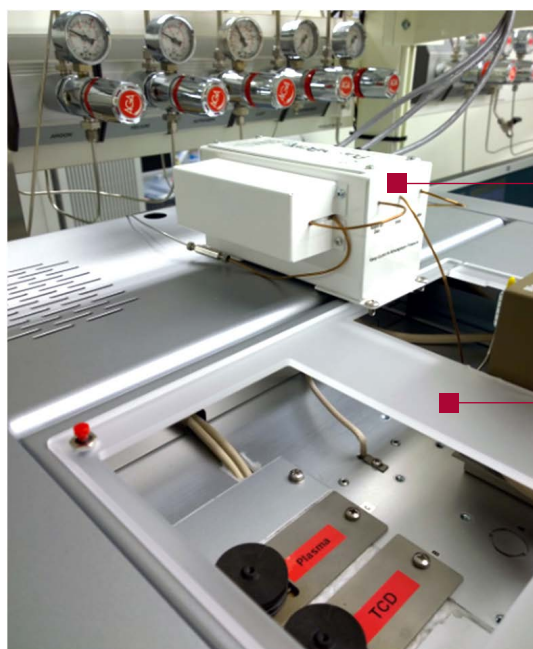

PlasmaDetek 2 (PED)



DANI Master GC

Koalent AB, Fältspatsgatan 6, 421 30 Västra Frölunda. Tel. +46 (0)31 709 81 80, e-post: koalent@koalent.se





CREATVAR
TECHNOLOGY CORP.

ANALYSIS OF TRACE IMPURITIES IN ELECTRONIC GRADE HYDROGEN USING AGILENT GC WITH PDHID AND PLASMADETEK 2

APPLICATION

Impurities(Ar-O₂-N₂-CH₄-CO-CO₂) analysis of electronic grade Hydrogen at 10 to 20ppb level with the PDHID and Plasmadetek 2.

This analyzer has been configured with two independent channels, front channel for impurities(Ar-O₂-N₂-CH₄-CO) measurement by PDHID and back channel for impurities(CO₂,n-C₂~4) measurement by PED. The system is based on a GC 6890N from Agilent equipped with electronic pressure control module. Two low leak rate 6 port and 4 port switching valves with Molsieve 5A column to perform bulk Hydrogen heart-cut technic and assure the baseline separation for Ar and O₂ in H₂ gas sample.

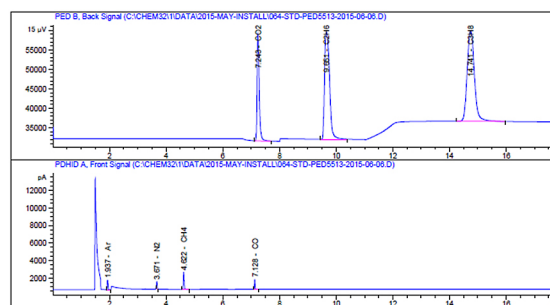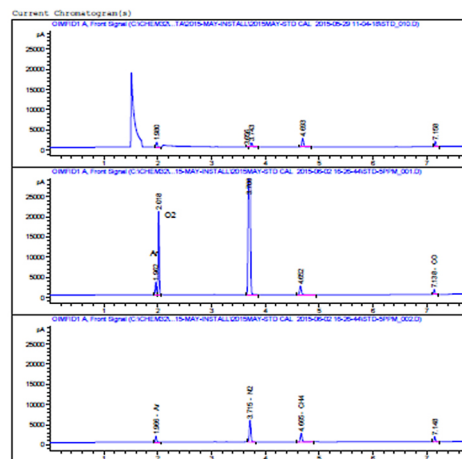

Signal 1: PED B, Back Signal

RetTime [min]	Type	Area [15 μ V*s]	Amt/Area	Amount [ppm]	Grp	Name
7.243	BB S	1.44087e5	2.78998e-5	4.02000	CO2	
9.651	BB S	3.36276e5	3.28599e-5	11.05000	C2H6	
14.741	BB	4.07393e5	1.13159e-5	4.61000	C3H8	

Signal 2: PDHID A, Front Signal

RetTime [min]	Type	Area [pA*s]	Amt/Area	Amount [ppm]	Grp	Name
1.937	BBA	1860.75464	2.86443e-3	5.33000	Ar	
3.671	BB	1753.38501	3.59305e-3	6.30000	N2	
4.622	BB	4349.97461	1.13334e-3	4.93000	CH4	
7.128	BB	1978.74536	2.41567e-3	4.78000	CO	

LDL CALCULATION				
COMPONENT	Standard Gas Conc. (ppm)	Peak Height	3xNoise	MDL (ppb)
AR-PDDA	5.33	1826.99624	0.3381	9.36858
N2-PDDA	6.3	794.99384	0.3381	2.679304
CH4-PDDA	4.93	1911.48582	0.3381	8.72009
CO-PDDA	4.78	1107.18896	0.3381	1.45966
CO2-PEDB	4.02	19544.75508	6.4824	1.33331
C2H6-PEDB	11.05	16159.39258	6.4824	4.43275
C3H8-PEDB	4.61	15624.92695	6.4824	1.91258



Copyright(C) 2008 CREATVAR TECHNOLOGY ALL Right Reserved
地址:807 高雄市三民區北平二街10號7F
Tel:07-3150125 Fax:07-3150025
www.cvtweb.com.tw



Agilent Technologies

Premier Solution Partner

Where innovation leads to success

www.ldetek.com

THE SMALLER, THE BETTER

LDP1000 SERIES



The new Compact-LDP1000 offers all features of the LDP1000, but in a smaller package:

- **2 BEDS OF PURIFICATION**

Allows perfect purification

- **RS-232 PORT**

Monitor the temperature of the 2 beds of purification

- **LEDS INDICATION**

Self-diagnostic and status of the purifier



When space is limited, it is the ideal solution for your purification need for noble gases, nitrogen and hydrogen. It is also possible to install it on the back panel of the MD2 to get a complete GC system in a single chassis.

Both LDP1000 and Compact-LDP1000 can give real end of life combined with the MD2 and PlasmaDetek technology. It is now easy to know when the getter needs to be replaced and avoid any problem with the system connected with the purifier.

AN ENHANCED DESIGN WITH **NEW CAPABILITIES**

MULTIDETEK 2



The MultiDetek 2 (MD2) has been well known in the industrial world with its rack mount design and plug and play philosophy, as well as being a unique compact GC design allowing multiple measurements and multiple components in one single chassis with easy maintenance.

By keeping this successful ideology, some features have been added to the MD2 to meet even more the needs of the laboratory world and achieve more applications. Compactness, versatility and flexibility is now possible for a single GC.

FOUR INDEPENDENT HEATED ZONES FOR VALVES/DETECTORS/SPLIT/SPLITLESS INJECTOR

The enhanced MD2 has now 4 independent heated zones for the valves and detectors. This design allows no cold point between each zone having a temperature up to 200 Celsius. The heated module covers are easy to remove and to put in place. Each module can contain up to 4 diaphragm valves with injector or one detector.

Where **innovation** leads to **success**

www.ldetek.com

1

CHROMALYTIC +61(0)3 9762 2034
ECHnology Pty Ltd

Website NEW : www.chromalytic.net.au E-mail : info@chromtech.net.au Tel: 03 9762 2034 . . . in AUSTRALIA

Australian Distributors
Importers & Manufacturers
www.chromtech.net.au

12 (of 16) 2016Mar

TWO INDEPENDENT LARGE RAMPING OVENS

Up to two larger ramping ovens can be installed on each side of the MD2 providing capabilities to use wider columns size as well as up to 16 columns in one oven. Temperature is also independent and can go up to 350 Celsius. An extra heated zone can be added on top of every large ramping oven for adding additional combination of valves and detector.



DataApex Clarity chromatography software compatible

COMPACT-LDP1000 MOUNTED ON THE BACK OF THE MD2

With the new design of the Compact-LDP1000, it is possible to mount the purifier on the back of the MD2 offering a complete solution on a single chassis. Isolation valve can be provided to isolate the purifier when needed.

HYDROGENICS

SHIFT POWER | ENERGIZE YOUR WORLD


Hydrogenics is a leading developer and manufacturer of hydrogen generation and hydrogen-based power systems. With a proven track-record of over 60-years of experience in fuel cell technology and innovation.

Hydrogenics helps companies reduce their footprint by reducing carbon emissions and producing clean renewable energy that can be sent back into the power grid. Hydrogenics super-dense PEM fuel cell allows for large-scale energy storage with unprecedented high-overload capability.

SOLUTION

The 1 MW pilot plant in Daesan Korea produces renewable electricity from unclean hydrogen and sends it back to the Korean power grid to be used for home and industry. Hydrogenics needed a company that offered advanced quality critical monitoring system and LDetek delivered with its turnkey solution cabinet system. The MultiDetek 2 checks the level of purity to ensure no contaminants in the hydrogen stream being delivered to the PEM fuel cell. LDetek provided a full system including cabinet, LDGSS (streams selector) and carrier gas purifier. As a company that values safety and reliability, Hydrogenics only chooses best-in-class partners, making LDetek an ideal choice to do business with.

INSTALLATION PICTURES

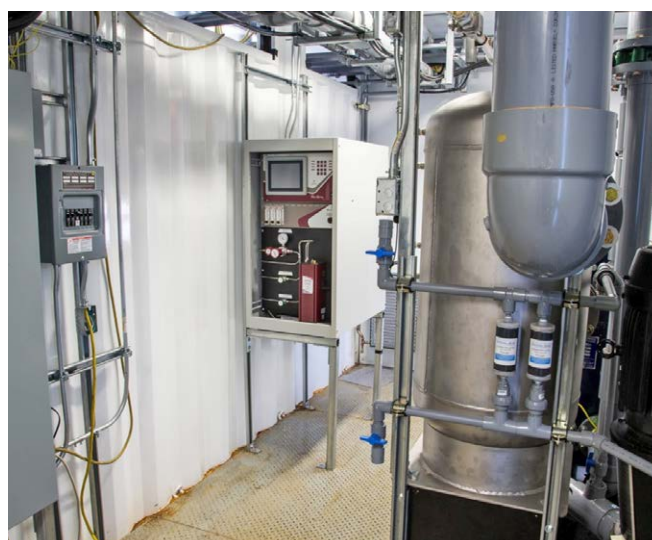

WWW.HYDROGENICS.COM

Where innovation leads to success

www.ldetek.com



COMPLETE PRODUCT LINE



MULTIDETEK SERIES

- **Compact Gas Chromatograph**
- One chassis configuration (6U Rackmount)
- Multichannels, Multimethods, Multidetectors
- Up to 6 isothermal or 3 programmable oven combination
- Up to 5 high purity proportional diaphragm valves (carrier-sample)
- Easy maintenance with its slide out design and front opening door
- ppb, ppm and % gas analysis
- Built in PC with 8.4" touch screen LCD & user-friendly interface
- Up to 10 high performance diaphragm valves
- Ethernet connectivity for remote control
- Integrated compact purifier with real end of life monitoring
- Serial/Profibus/Modbus communication protocols
- Fast parallel chromatography
- Multi heated zones to avoid cold points
- Purged & real time monitored zones for hazardous gases
- Multi sample injection techniques



LABDETEK

- **Stand alone heated zone module**
- Large fast programmable oven up to 350 Celsius
- Isothermal zone up to 200 Celsius
- Multiple columns, valves, detectors, I/Os
- Integrated gas purifier system
- Multiple flow controllers
- User-friendly interface to control the device.
- Compatible with Clarity from Data Apex



PLASMADETEK SERIES

- **Stand alone detector for any GC**
- Heating up to 200 Celsius
- Ideal for ppb-ppm trace impurities
- Selective, sensitive and configurable
- Replaces ECD-FID-TCD-DID all in one detector
- Compatible with Helium, Argon, Nitrogen and Hydrogen as carrier gas
- Interchangeable optics & carrier gas type within seconds
- Compatible with Clarity from Data Apex



LDP1000 SERIES

- **Purification of Noble gases, Nitrogen and Hydrogen**
- Interchangeable getter philosophy
- Enhanced 2 beds of purification
- Different gas flow capability
- Real end of life monitoring in combination with the PlasmaDetek series and MultiDetek 2



LD8000 SERIES

- **ppb/ppm N2 in Ar - He - Ne - crude argon - multigas**
- Built-in multi-stream selector system.
- Integrated zero calibration gas generator system.



LD2000 SERIES

- **ppb/ppm THC in air - O2 - N2 - CO - CO2 - H2 He - Ar - Ne - Kr - Xe**
- Electronic flow controllers for air, fuel and sample.
- Safety fuel shutoff valve.



LDGSS SERIES

- **Stream selector system for vacuum up to high pressure gas lines**
- Options available for corrosive and toxic gases
- Leak proof, no stream contamination, no outboard leakage, no cross contamination
- Can be connected to the MD series for remote control switching and sequence programming
- Double block and bleed version available



Be sure to download our
CHROMATOMAG
 and visit our application notes section at

www.ldetek.com



T