




### AUTOSAMPLER QUEUE and BATCH REPROCESSING

Create customized autosampler sequences for liquid injections, purge and trap autosamplers, gas sampling valves, and stream selectors including unique predefined sample information, auto-calibration and batch reprocessing of previously run samples.

### DATA LOGGER MODE

Peak Simple's Data Logger Mode allows you to display a scaled and calibrated result in large numbers instead of the usual strip chart data presentation. Data Logger Mode is especially useful when monitoring total hydrocarbons on one channel while performing a separation on another channel.

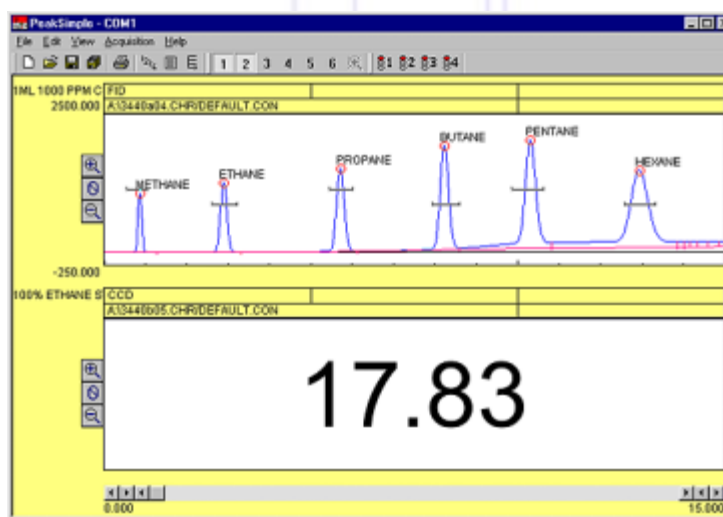

### SELF-VALIDATING HARDWARE

PeakSimple will play back and re-acquire any chromatogram multiple times, establishing the precision and accuracy of the data system using real data, not "canned" chromatograms. PeakSimple's validation can be performed by the user anytime without extra hardware.

