



Bottle not included.

#### Waste Overflow Indicator for LC Systems

- Avoid messy pooling around mobile phase waste containers.
- Audible alarm instantly alerts user, preventing overflow.
- Compact, battery-operated unit.
- Available for 4-liter and GL-45 solvent bottles.

The Restek® waste overflow indicator will help keep your mobile phase waste where it belongs—in the waste container! Compact, battery-operated unit fits securely on solvent bottles and accommodates two waste streams. An audible alarm is given as the solvent waste container approaches capacity, giving you time to empty or change the container. Another innovative design from Restek!

Description	qty.	cat.#
Waste Overflow Indicator for LC Systems, 4 Liter	ea.	26543
Waste Overflow Indicator for LC Systems, GL-45	ea.	26550
Replacement AA Battery for the Waste Overflow Indicator	ea.	26544
Replacement AA Batteries for the Waste Overflow Indicator	3-pk.	26545

Extend the life of columns and pump seals with Restek's Bluestem glass solvent filter

See **page 343**.



#### tech tip

#### Mobile Phase Additives with LC-MS

Mobile phase additives such as triethylamine, trifluoroacetic acid, and ion-pairing reagents can compete with sample ions, decreasing sensitivity and, in some cases, reducing sample ion intake into the MS. To obtain symmetric peaks and/or sufficient retention, use base deactivated silica, like Pinnacle® DB or high purity silicas like Ultra or Raptor™, that minimize the need for additives.



26395



Polypropylene Membrane Filters



25925

#### Hub-Cap Filter Kit


Kit includes: bottle adaptor, bottle adaptor nut, filter inlet cap, grid support, vacuum hose barb, tube compression fitting, 47 mm grid, 47 mm 0.22 µm filter membrane, 47 mm 0.45 µm filter membrane, 1/4" OD x 1/8" ID ultra chemical resistant, FEP-lined Tygon® tubing (3'), 6" x 6" box with shrink-wrap insert.

Includes universal threads designed for 4 L or Wheaton bottles. Use with GL-45 bottles requires Opti-Cap® adaptor (sold separately).



Bottle not included.

#### Assembles quickly and easily!

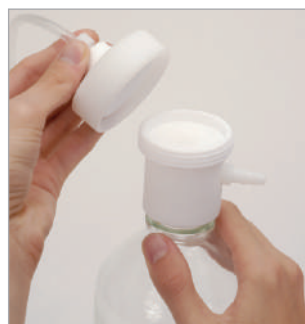

Unscrew and lift off top.



Place membrane filter on top of grid.



Reattach top.



Connect vacuum line to side port.

Description	qty.	cat.#
Hub-Cap Filter Kit for 4 L or Wheaton bottles	kit	26395
Replacement Parts	qty.	cat.#
Polypropylene Membrane Filters, 47 mm, 0.45 µm	100-pk.	26396
Polypropylene Membrane Filters, 47 mm, 0.22 µm	100-pk.	26397
Nylon Membrane Filters, 47 mm, 0.45 µm	100-pk.	26398
Nylon Membrane Filters, 47 mm, 0.22 µm	100-pk.	26399
Hub-Cap Filter Hose Barb	ea.	25925